

Multi-Adjustable



“The adaptable chair solution,
designed for multiple users”



Dr Hilary Jones
Our In-house Health Advisor

COMFORT WITHOUT COMPROMISE



**Considered users,
those who require:**

- Postural support
- Pressure care management
- Assistance from sitting to a stand position
- High leg elevation

Suitable for:



Home



Healthcare



Nursing and
Residential Homes

The Multi chair range is an innovation in dynamic frame and chair design providing a cost effective solution to the demands of changing care needs, infection control and efficiency.

These chairs can be fully adapted by the healthcare team in the home or on the ward. Each chair can be quickly and easily customised for different users or for an individual as their care needs change over time.

The Multi frame and full range of compatible cushions gives you complete flexibility and control. You choose the right seat height and seat cushion, the right back rest and lateral support for the individual you are caring for. Adapt the chair over time to provide more support and pressure care or redesign the chair for a new person. The Multi range just keeps on working.

The Multi C-air is our standard Multi rise and recline chair we also have Bariatric and Porter styles in the Multi range.



COMFORT WITHOUT COMPROMISE

Multi-C-air

The Multi C-air an innovative and dynamic design that solves the challenge of continual changing care demands.

- Adapt the Multi C-air for any size or shape of individual with interchangeable waterfall or lateral back cushions.
- Support different posture requirements and positions with alternative and additional lower and upper back and neck cushions.
- Manage specific comfort and pressure management needs with interchangeable seat cushions constructed from a range of specialist pressure care materials.
- Control infections by choosing from our range of healthcare grade fabrics.
- Maintain the highest hygiene standards with quick change covers and replacement cushions.
- Purchase replacement or new support cushions as you need them to extend the use of your Multi C-air.

Key features

- 1 Separate waterfall back cushions can be individually repositioned for comfort and support.
- 2 Back cushions can be removed and replaced with more supportive lateral styles.
- 3 Reduce the seat depth with adjuster pads that fit behind the back cushions.
- 4 Choose internal seat cushions from a range of pressure care management solutions.
- 5 The seat and leg chaise are generously padded with memory foam for comfort and support.
- 6 To fit the Multi C-air to a taller person you can use the 2" leg adjusters to raise the seat height and adjust the chair with a longer leg rest.
- 7 Back and seat cushions are upholstered in a vapour permeable multi-way stretch fabric, to aid pressure care. The fabric is washable at 95°C and available in a range of colours.
- 8 Glides are built-in as standard to the Multi C-air.
- 9 Single or dual tilt in space mechanisms are available for the Multi C-air.



Leg Height Adjusters




The **Multi C-air**, the hardest working rise and recline chair available.
An **effective, efficient** and **value for money** solution.

Repose Multi C-air seat options for pressure management

All cushions are covered in multi-way stretch, vapour permeable fabric which is wipeable, removable, washable at 95°C, incontinence proof and breathable.

Viscoform

- Memory foam (std)
- Reflex foam

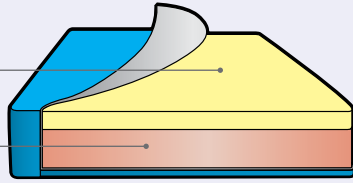

Pressure Risk


LOW MEDIUM

For those who have good posture and mobility but are seated for greater than one hour at a time and require pressure relief and supported comfort.

This cushion has a high density foam base which offers a supportive and firm seat base with a top layer of memory foam to allow submersion into the cushion giving stability and a large surface area to reduce pressure levels.

Liquiform

- Memory foam
- Liquid gel sac
- Reflex foam

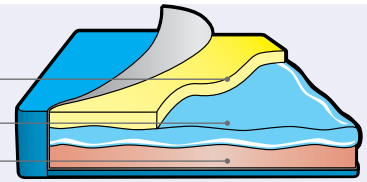

Pressure Risk


MEDIUM HIGH

For those who have good posture and mobility but are seated for long periods of time and require pressure relief and supported comfort. Excellent for people with little fatty tissue.

This cushion has a high density foam base which offers a supportive and firm seat base topped with a fluid liquid gel in a soft silicone cover to provide high levels of pressure relief. The liquid gel helps to dissipate heat and constantly adapts to movement.

Celliform

- Gelacell mesh gel
- Reflex foam

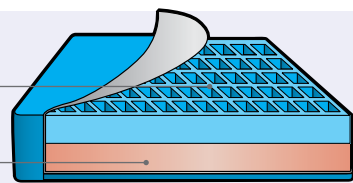

Pressure Risk


MEDIUM HIGH

For those sitting for long periods of time with specific comfort issues such as a high risk of pressure sores and shear forces.

Gelacell's unique grid design eliminates pressure points by contouring in response to the transfer of weight and evenly distributing pressure. The molecular structure restricts absorbed heat and keeps the seat cooler and fresher.

Airform

- Memory foam
- 3 x manual air cushions
- Reflex foam

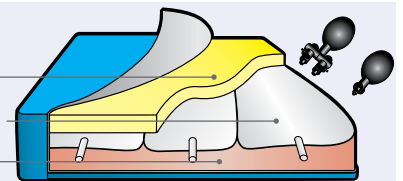

Pressure Risk


HIGH

For those who have good postural position but who are unable to alter their own position easily and require high level pressure care.

This cushion gives the benefit of an air cushion without the sensation of alternating changes in air pressure.

Dynaform

- Westmeria alternating air cushion
- Reflex foam

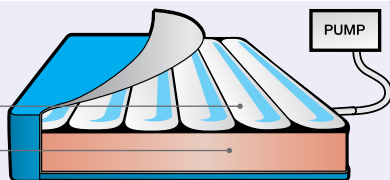

Pressure Risk

VERY HIGH

For those at very high risk of developing pressure sores or a history of sores where the highest priority is to prevent further pressure damage.

For users only up to 20 stones (120kg). The Dynaform cushion is a dynamic alternating air cushion which constantly changes the air levels within its many air cells to recreate the effect of pressure relieving and prevent pressure build up.

Air cell features


- 6 dynamic cells 4" (10cm) wide
- 3 cell cycle repeated every 9 minutes

Please note the Dynaform cushion does need an electric supply at all times to operate.


Pressure mapping

All Repose seat designs undergo rigorous pressure mapping testing to ensure we achieve the highest standards possible. A mat of pressure sensors is applied to the seat. The matrix of sensors quantifies the surface pressure between a person and the support surface the chair provides. The results provide a map for analysis by clinicians who can then advise our design teams. The aim is to ensure that pressure is as evenly distributed as possible and that the maximum area of the thighs and buttocks is supported as at low an interface pressure as possible.

Poor pressure distribution

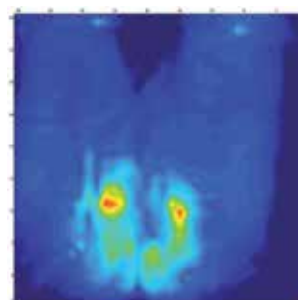

Good pressure distribution

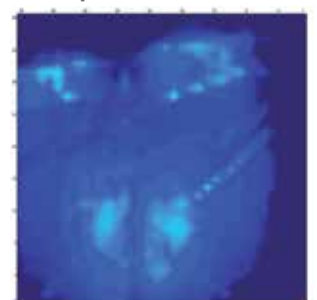


Pressure Map for individual: Gender: Female; Age: 35 Years; Weight: 10.5 Stones.

Note: Users with pressure relief needs or concerns are advised to seek the advice of a tissue viability clinician on the most suitable level of pressure relief for their medical requirements.


COMFORT WITHOUT COMPROMISE

Back cushion options for enhanced support

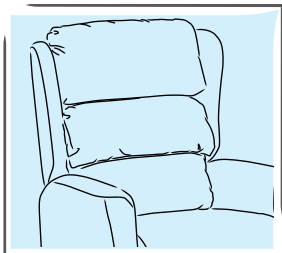


Waterfall back cushion

Individual waterfall cushions, which are fully adjustable to accommodate maximum postural support.

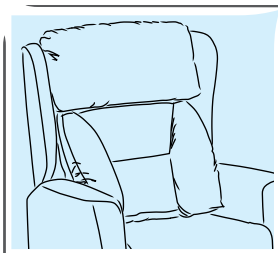


Comfort lateral back cushion

Soft fibre filled removable lateral cushion for lower back and side lumbar support, without feeling restricted.

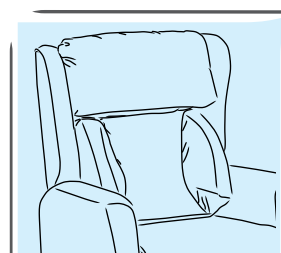


Support lateral back cushion

A foam filled back providing firmer lateral support with a soft fibre filled headrest.

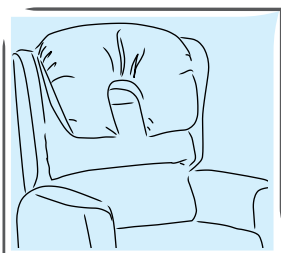


Large profile headrest with waterfall back

Recommended for those with limited head control and requiring a soft adjustable lumbar support.

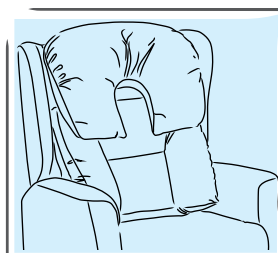

Large profile headrest with comfort lateral back

Recommended for those with limited head control and requiring soft lateral support without feeling restricted.

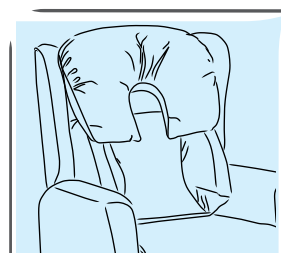

Large profile headrest with support lateral back

Recommended for those with limited head control and requiring firm foam lateral lumbar and midline support.

Care features



Easy replace covers. Quick to change and replace for care and cleaning.



Interchangeable seat cushions to adapt the Multi C-air for each new user's pressure care or support needs.



Flexible positioning of the three waterfall cushions to get the shape of the back right for each user, with adjustable back padding in all cushions.

Optional extras

- Head roll headrest
- Small profile headrest.
- Standard arm caps and full arm caps.
- Honeycomb gel in arm rests.
- Replacement seat chaise cover and memory foam.

More options are available, contact for more information

Rise and recline mechanisms

The Multi C-air comes with the following rise and recline options:

Single Motor – Tilt in Space




Full tilt




Static position



Semi recline

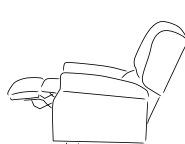


Footrest horizontal

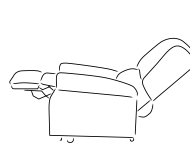

Fully reclined


Dual Motor – Tilt in Space

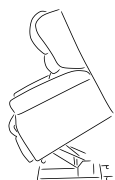


Full tilt




Static position

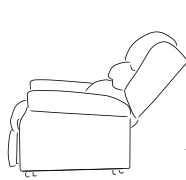


Independent back rest

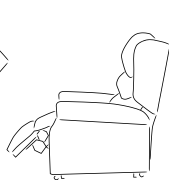

Semi recline

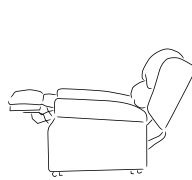

Footrest horizontal

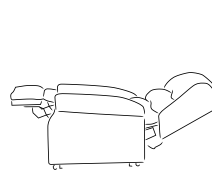

Fully reclined

Available sizes

	SMALL	MEDIUM	LARGE
SEAT HEIGHT	16" (41cm)	18" (46cm)	20" (51cm)
SEAT WIDTH	18" (46cm)	20" (51cm)	22" (56cm)
SEAT DEPTH	18" (46cm)	20" (51cm)	20" (51cm)
By fitting 2" (5cm) leg adjusters the SEAT HEIGHT is 2" higher	18" (46cm)	20" (51cm)	22" (56cm)
By fitting 2" (5cm) seat depth pad the SEAT DEPTH is reduced by 2"	16" (41cm)	18" (46cm)	18" (46cm)

Weight limit

The Multi C-air has a weight limit of 20 stones as standard. A 25 stone option is available please ask for details.

Fabric options

- Waterproof
- Anti MRSA
- Anti Bacterial
- Anti Fungal
- Crib 5
- Bleach cleanable
- Vinyl.

The fabric options shown in this brochure: SMD Bolsena Dewberry with a Plain Blue vapour permeable stretch seat and back cushion covers.



Guarantees

The Repose Multi C-air comes with the following guarantees:

- 2 Year electrical warranty for the recliner mechanism.
- 1 Year warranty on alternating seat cushions.
- 1 Year warranty on the upholstery and frame.*

If you have any questions about the Multi C-air or the Multi range, please contact us on: 0844 7766001 or at info@reposefurniture.co.uk

*Full terms and conditions available upon request



Thank you for supporting British manufacturers

Repose Furniture Ltd reserves the right to modify or change the specification without prior notification. For terms and conditions of sale and or agreement please refer to your supplier.



Specialist
Seating Solutions
Repose Furniture Ltd.
Attwood House, Cokeland Place
Cradley Heath, West Midlands
B64 6AN
Tel: 0844 7766001
Email: info@reposefurniture.co.uk

Recommended supplier: